

AZURA® pH flow cell

The AZURA® pH flow cell is designed for both micro and semi-preparative applications.

Combined with an inserted AZURA® pH electrode, which is connected to the AZURA® Conductivity Monitor CM 2.1S, it allows continuous pH monitoring of fluid streams in FPLC systems or other flow applications.

The AZURA® pH flow cell is delivered with various mounting screws and capillary-fitting sets to allow flexible assignment into an AZURA® FPLC system according to your needs.

Scope of delivery

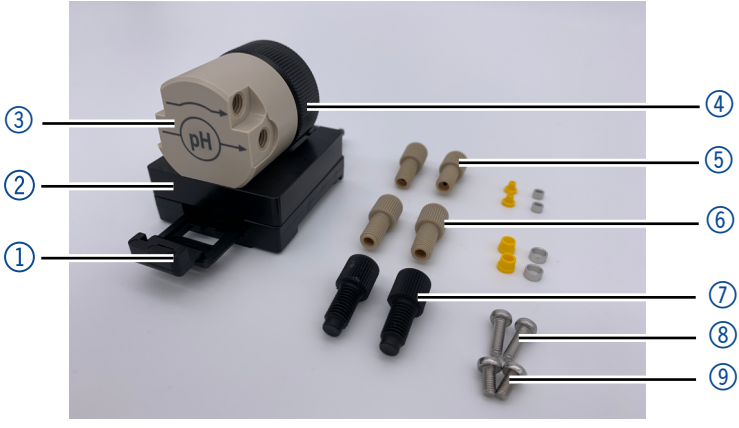
Legend	
① Universal AZURA® Click mounting rail holder	 <p>The image shows the AZURA® pH flow cell assembly. On the left, the flow cell (3) is mounted on a black adapter plate (2) which is attached to a universal click mounting rail holder (1). The flow cell has a 'pH' label and a cap nut (4). To the right, various components are laid out: a capillary fitting set for 1/16" tubing (5), a capillary fitting set for 1/8" tubing (6), two black blind fittings (7), two screws of size M3x6 (8), and two screws of size M3x10 (9).</p>
② Adapter plate	
③ AZURA® pH flow cell	
④ Cap nut of the AZURA® pH flow cell	
⑤ Capillary fitting set for 1/16" tubing	
⑥ Capillary fitting set for 1/8" tubing	
⑦ Blind fittings	
⑧ Screws, size M3x6	
⑨ Screws, size M3x10	



Fig. 1: AZURA® pH flow cell, scope of delivery



Mounting


The AZURA® pH flow cell ③ can be implemented into the AZURA® FPLC system via the AZURA® Click system, the AZURA® Organizer ⑩ or the AZURA Mounting Bracket Bio L ⑪ . For further information on how to install the flow cell refer to the supplementary sheet of the assembly of your choice (www.knauer.net/library).

Note: In order to mount the AZURA pH flow cell onto the AZURA® Organizer or Mounting Bracket Bio L, the universal AZURA® Click mounting rail holder ① has to be unscrewed first.

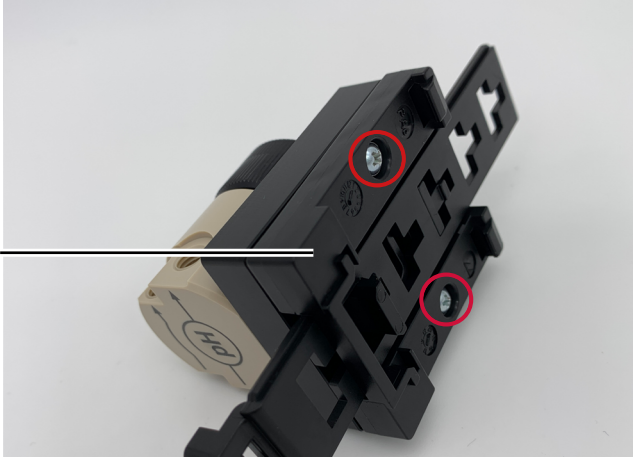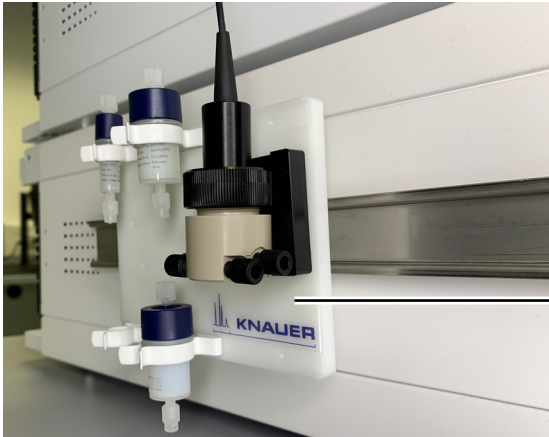- Tools**
- Phillips screwdriver, PH1
 - Torx screwdriver, T10
 - 2 screws, size M3x06
 - 2 screws, size M3x10

Process	Figure
<ol style="list-style-type: none"> 1. Loosen the two screws on the back of the assembly with a Phillips screwdriver (PH1). 2. Remove the universal AZURA® Click mounting rail holder ① . 	 <p>Fig. 2: Loosen two screws</p>
<ol style="list-style-type: none"> 3. Now you can mount the pH flow cell to the AZURA® Organizer ⑩ by using the two screws (size M3x10) and the T10 torx screwdriver. 	 <p>Fig. 3: pH flow cell mounted to AZURA® Organizer</p>

Process	Figure
<p>4. Or you can mount the pH flow cell to the AZURA® Mounting bracket Bio L ⑪ by using the two screws (size M3x06) and the T10 torx screwdriver.</p>	 <p>Fig. 4: pH flow cell mounted to AZURA® Mounting bracket Bio L</p>
<p>5. If you want to mount the pH flow cell to the AZURA® Click mounting rail ⑫, you do not have to remove the AZURA® Click mounting rail holder.</p>	 <p>Fig. 5: pH flow cell mounted to AZURA® Click system</p>

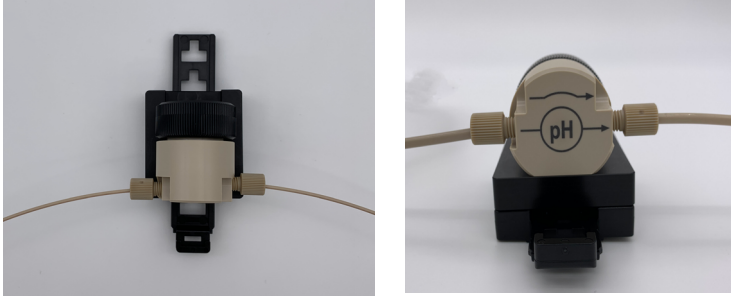Process	Figure
<p>6. Connect the AZURA® pH flow cell to the liquid capillary by using the included flat-bottom fitting set with either 1/16 or 1/8" OD polymer capillaries.</p>	 <p>Fig. 6: AZURA® pH flow cell with connected liquid capillaries, (1/16" or 1/8" OD)</p>

Operation

Inserting a pH electrode

To implement an AZURA® pH electrode, unlock the cap nut ④ from the pH flow cell ③. Insert the pH electrode and screw the cap nut hand tight back onto the pH flow cell.

Next steps:

Connect the BNC plug to the BNC connector of the CM 2.1S.

For further information refer to the instructions of the AZURA® Conductivity Monitor CM 2.1S.



Note: Always ensure not to build up a pressure inside the pH flow cell that exceeds 5 bar. The max. pressure rating of the AZURA® pH flow cell is higher than the burst pressure of classic pH electrodes, which is usually around 5 bar. The max. pressure of the bypass is only limited by the utilized capillary-fitting set.

Removing a pH electrode

In order to remove the pH electrode, unlock the cap nut and store the pH electrode in a buffered solution as recommended in the corresponding instruction manual.

Inserting and removing a dummy electrode

Lock or unlock, respectively the cap nut. Clean the dummy electrode with water and store dry in a clean environment.

Tips for FPLC system operation without inserted pH electrode

pH electrodes have to be removed from the system flow, either because they:

- have to be replaced or recalibrated
- do not withstand the used eluent or cleaning solvent

In order to still be able to operate the FPLC system you can either:

- connect the fluid capillaries to the integrated bypass (use blind-fittings during the tubing switching to prevent leakages, see figure 7) or
- insert the dummy electrode (order separately) after removing the pH electrode (see figure 8).

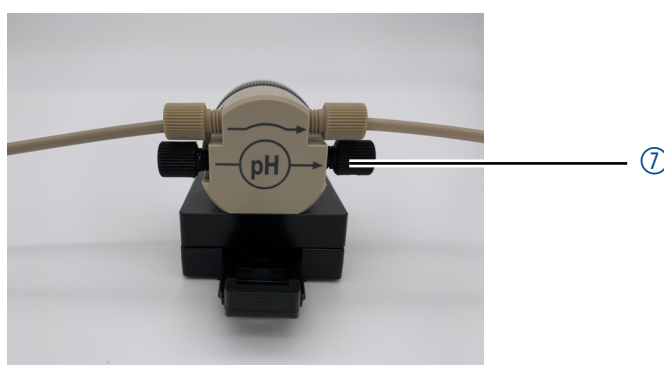


Fig. 7: AZURA® pH flow cell with blind fittings

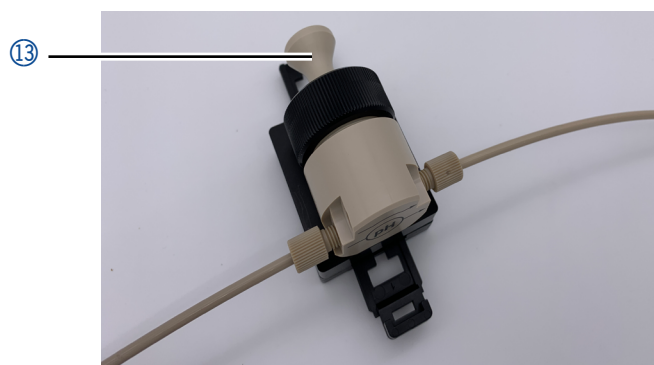


Fig. 8: AZURA® pH flow cell with dummy electrode

When utilizing the integrated bypass, the overall system volume minimizes. To keep the dwell volume constant (for fractionation purposes), utilize the AZURA® pH dummy electrode ⑬ instead.

Repeat orders

Name	Order Number
pH electrode, round tip, with quality control certificate	A1933-1
AZURA® pH flow cell incl. mounting set, 1/4-28" UNF FB fittings for 1/16" & 1/8" OD tubing and blind fittings	A1943
AZURA® pH dummy electrode, PEEK, incl. cap for AZURA pH sensor (A70091-2)	A1942-1
O-ring for A1943, FFKM	A1943-1